Taiwan
Advanced
Nanotech

(Compatible with power 220V-240V only)

Maelstrom 9610

User Manual

TAN Bead Nucleic Acid Extractor
(Non-sterile)

039.M1601.X31

Table of Contents

1.	Introduction	6
2.	Instrument Overview	13
3.	Installation & Get Started	17
4.	User Interface Guidelines.....	18
5.	Technical Support	40
6.	Cleaning and Maintenance	41
7.	Disposal	41
8.	Patents	42
9.	About Manufacturer	42
10.	Appendix	43

About Manual

The label on the instrument, the User Manual, and other packaging material may contain following symbols:

Catalog number

Non-Sterile

Serial number

Temp Limit
(Operating)

Manufacturer

Humidity Limit
(Operating)

The date of manufacture

Biologic Risks

Keep Away From Sunlight

Keep Dry

Consult Instructions for Use

The CE mark

In Vitro Diagnostic Medical Device

Authorized European Representative

Caution, consult accompanying documents

Hot surface, risk of burns

Watch your fingers and your hands

The symbol indicates that the product should not be discarded as unsorted waste but must sent to separate collection facilities for recovery and recycling

About Instrument

WARNING

- Use a power cord that meets your country's standard. In case of any questions about the power cord, contact your local distributor or reseller.
- Do not use the instrument with voltage other than the voltage specified on the device.
- Do not use the instrument with damaged power cord or a loose socket.
- When the power plug is disconnected from the AC outlet, please hold the power plug itself instead of pulling the power cord.
- Before working on maintenance, must disconnect the power plug from the outlet in advance.
- Do not pour any liquid on the instrument.
- Do not place any containers with liquid on the instrument. Doing so may cause a fire, electric shock or malfunctions of the instrument.
- Do not touch the power plug or cord if there is a chance of lightning. Failure to observe this may cause electric shocks.
- If you hear thunder or suspect approaching lightning, quickly turn off the power switch and pull the power plug from the AC outlet. Failure to observe this may cause a fire or malfunctions.

About Instrument

CAUTION

- Never attempt to remodel the instrument without the permission from the manufacturer. Doing so may lead to a fire or electric shock.
- Do not bump or knock the instrument. Doing so may cause malfunctions.
- Repairs to the instrument must be performed by agencies authorized by Taiwan Advanced Nanotech Inc.
- Only original spare parts from Taiwan Advanced Nanotech Inc can be used on the instrument.
- If the equipment is used in a manner not specified by the manufacturer, the protection given by the instrument may be impaired or invalid.


1. Introduction

Maelstrom 9610 is an automated nucleic acid platform designed for high-throughput applications. Specialized spin tips enable super mixing efficiency of magnetic beads and larger processing volume. With intuitive interface and flexible program, Maelstrom 9610 can boost molecular diagnostic laboratory productivity by transforming routine operations into a walk-away solution.

Operation Principle

Maelstrom 9610 uses magnetic rods to collect and transfer magnetic beads and spin tips to mix suspension. Rotating the spin tips will mix suspension very well. Extended magnetic rods can collect magnetic beads efficiently.

Purified nucleic acids can be obtained after cell lysis, nucleic acids absorption, wash and elution.

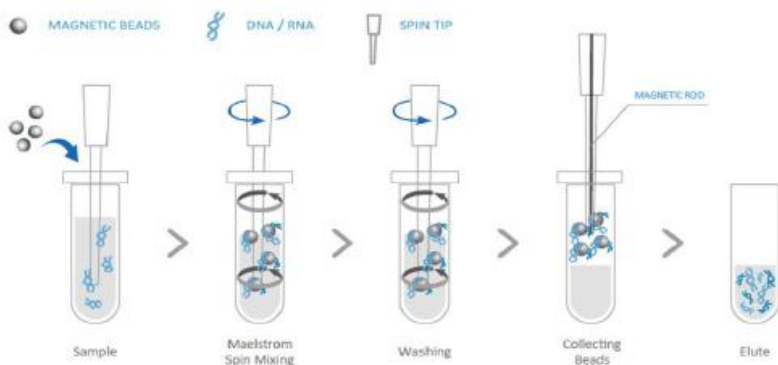- The gear used in the TANBead Nucleic Acid Extractor is a hollow gear, which allows the magnetic rod to pass through the connecting portion to insert into the spin tips to perform the magnetic effect. In addition, the gear comprises an extended hollow bearing, such that the collimation of the stirring device is increased.

1. Introduction

- The stirring device performs mixing by using the characteristics of rotation, biggest flux per unit area can be achieved compared to up and down.
- The rotation stirring design can decrease the air mixed into the sample and avoid problem of cross contamination in experimental.

1. Introduction

Intended Purpose

Maelstrom 9610 automated nucleic acid extractor is an in vitro diagnostic device, which is used to extract and purified nucleic acids from human and biological samples. The purified nucleic acid can be used with any downstream application employing PCR-based qualitative, semi-quantitative and quantitative assays. This product should be operated by technicians, physicians, and biologists with well-trained in molecular biological techniques. It is suitable for all population in human.

Maelstrom 9610 automated nucleic acid extractor is equipped with different types of TANBead Nucleic Acid Extraction kits, which can extract and purify nucleic acids from biological samples such as blood, cfDNA, virus, tissues, Bacteria and stool. Please refer to the following for types of specimens required.

Sample Types

Maelstrom 9610 Nucleic Acid Extractor is used with TANBead® nucleic acid extraction reagent kit. It would be suitable for the model names as below. Please refer to the instruction of TANBead Nucleic Acid Extraction kit.

Sample	Blood			cfDNA	
Kit	611 series	61E series	621 series	61C series	

Sample	Virus				
Kit	615 series	635 series	665 series	685 series	61H series

Sample	Tissue			Bacteria	Stool
Kit	612 series	6T2 series	6K2 series	61G series	6SC series

1. Introduction

Environmental Requirements

To avoid shortening of the lifespan of the instrument, use Maelstrom 9610 in a location that meets the following criteria:

- Choose a location that has good air circulation.
- Place Maelstrom 9610 on the table capable of bearing at least the weight of Maelstrom 9610 (≈ 95 kg).
- Do not use Maelstrom 9610 in a location where there is large temperature variability or high humidity.
- Suggested operation condition:
Temperature: 10-35°C
Relative humidity: 40-80%
- Storage and transportation condition:
Temperature: 5-60°C
Relative humidity: 20-90%
- Maximum operating altitude: 2000 m

1. Introduction

Safety Instructions and Guidelines

- This device can be used with potentially biohazardous materials. Use appropriate personal protective equipment (gloves, safety goggles, lab coat, etc.) for handling and disposing of biohazardous materials.
- Under normal condition, sound pressure level from Maelstrom 9610 does not exceed 80 dB and cause a hazard. Please contact technical support in case of higher sound pressure level.
- This device can be hazardous due to the use of chemical and biohazardous substances.
- Users should adhere to their institutional guidelines for the handling and disposal of all infectious substances used with this device.
- It is important to clean the device after every use. If samples or reagents have been spilled, it is important to clean the instrument immediately to avoid damage or contamination of samples.
- This device is to use with the compatible spin tips. Using incompatible spin tips may cause poor extraction performance.
- Read this user manual in its entirety prior to operating the device. Failure to read, understand, and follow the instructions in the manual may result in damage to the device, injury to laboratory and operating personnel or poor performance.

1. Introduction

Safety and EMC Information

Safety Requirements

The medical device has passed the tests and conformed to the following standards;

- IEC 61010-1:2010+A1:2016 (Edition 3.1) and/or EN 61010-1:2010+A1:2019 / DIN EN 61010-1:2020, "electrical equipment for measurement, control, and laboratory use Part 1: General requirements"
- IEC 61010-2-101:2018 with IEC 61010-1:2010 + A1:2016 / EN IEC 61010-2-101:2022+A11:2022 with IEC 61010-2-101:2018 / EN 61010-2-101:2017 / DIN EN IEC 61010-2-101:2023 "electrical equipment for measurement, control, and laboratory use Part 2-101: Particular requirements for in vitro diagnostic (IVD) medical equipment"

EMC Requirements

The medical device has passed the tests and conformed to the following standards;

- IEC 61326-1:2020 / EN IEC 61326-1:2021 "Electrical equipment for measurement, control and laboratory use – EMC requirements - Part 1: General requirements"
- IEC 61326-2-6:2020 / EN IEC 61326-2-6:2021 "Electrical equipment for measurement, control and laboratory use – EMC requirements - Part 2-6: Particular requirements - In vitro diagnostic (IVD) medical equipment"

1. Introduction

Equipment list

Please note that following accessories accompanying with instrument may vary from region to region.

Power cord x1

Transmission cable x1

Magnetic separator x1

Tips combo pack x1

Auto Tube base x 6

User manual x1

=Remark=

Power cord varies from region to region.

Auto Tube base is required in use of Auto Tube reagent kits.

2. Instrument Overview

ITEM	SPECIFICATION
Model	Maelstrom 9610
Weight	Approx. 95 kg
Dimensions	870(W) x 575(L) x 700(H) mm
Power rating	220-240Vac, 50/60Hz, 3.5A (Class I)
Fuse	250V, 5A
Max. Throughput	96 samples per run
Process. volume	50µl -1,600µl
Spin speed	up to 3000 rpm
Heater	4 independent heating plates
Magnetic rod	>3900 gauss
Display	7-inch touchscreen
UV & HEPA	UV: UV-C type 4W, HEPA: E10 Class
Program Name	Firmware: PLC; Adjustment Program: C#

Please note that this machine is compatible with 220-240V power only. Using wrong power may lead to malfunction or damage to the instrument.

2. Instrument Overview

USB type A


USB type A port is for the use of barcode scanner.

USB type B

USB type B port is for the connection with PC.


① Power Button & USB I/O

② Turntable rotating buttons

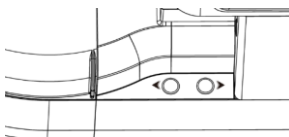Use "►" and "◄" button to rotate turntable.


③ Turntable

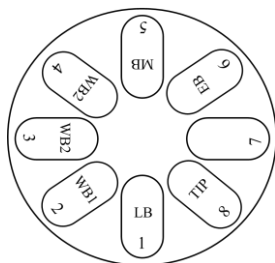- Turntable contains 8 plate holders.
- There are four plate holders (NO.1, NO.2, NO.3 and NO.6) have the heating module for temperature control.
- No. 8 plate holder is for mounting spin tips.
- NO.2 to NO.5 plate holder are allowed to be selected as ejecting spin tips.

2. Instrument Overview

Rear View

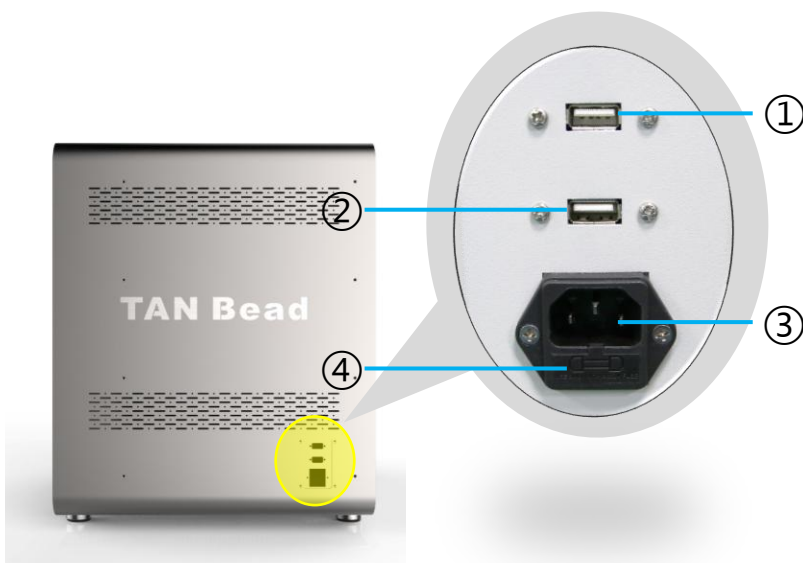- ① USB type A for use of upgrading HMI program
- ② USB type A for use of upgrading PLC program
- ③ Power inlet
- ④ Fuse

Please note that only TANBead technicians and authorized technicians or service engineers shall engage in HMI and PLC program upgrade.

2. Instrument Overview

Control Panel

To turn off the alarm

To pause the program. Press again to resume the program.

To turn on/off the LED light.

To initiate the program

Emergency to stop machine

3. Installation & Get Started

Installing Maelstrom 9610 is easy. Please follow the guide below to correctly set up the instrument.

Please note that this instrument weights around 95kg, it is highly recommended to have 3 or more persons to handle when taking it out from the box and must pay attention to the safety.

Step 1:

Take out the instrument from the box (3 or more persons to handle is highly advised)

Step 2:

Remove the cushions around the instrument

Step 3:

Place the instrument to a flat-table that can bear 120kg weight

Step 4:

Connect the power. Please note that this instrument is compatible with 220-240V power only. Using wrong power source will lead to malfunction or damage to the instrument

Step 5:

Power on the instrument and it will perform initialization. Once initialization is done, please remove Y-axis protector from the instrument and the instrument is ready to use.

◀ Y-axis protector

4. User Interface Guidelines

Create an User ID

An user ID is required to access the instrument. For a brand new machine, there is a need to create an user in advance. Click the icon to create your own user ID.

◀ Press the icon to create a new user ID

◀ Click blank space to call out keyboard

◀ Use the keyboard to input user ID and then click [Enter] to save and exit

User ID: max. 12 characters


4. User Interface Guidelines

Save User ID

After creating the user ID, please must click  to save and then program will automatically lead you access the system.

◀ Instrument is trying to save new user ID

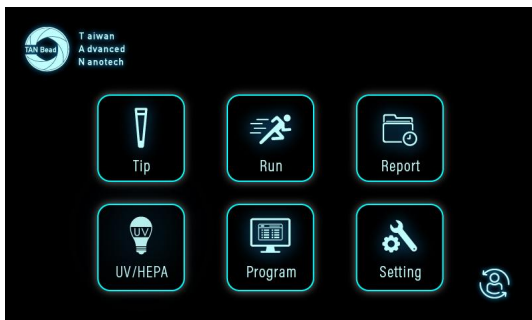◀ You will be guided to the main page after saving. Now you can prepare for extraction work.

◀ Switch user

To switch user, there is no need to power off/on the instrument. Take advantage of icon  to switch user quickly.

4. User Interface Guidelines

Possible Error Message

Some messages may occur when you try to log in or create an user ID.

- ◀ **User ID not found**
Make sure you input correct, existing user ID account

- ◀ **User ID repeated**
Please change the name of creating an ID

- ◀ **Account is full**
Max. 24 user ID are allowed. If the account is full, please contact the dealer nearest to you for further help.

4. User Interface Guidelines

Main Page

There are six functions that can be used on main page.

To test on mounting and ejecting spin tips

To execute the extraction program

To export extraction records

To use UV & HEPA to clean the instrument

To edit, add, and delete extraction program


To set basic parameter of instrument

4. User Interface Guidelines

Tip (Test on Mount and Eject Spin Tips)

Please complete the tip mounting before you start the extraction protocol.

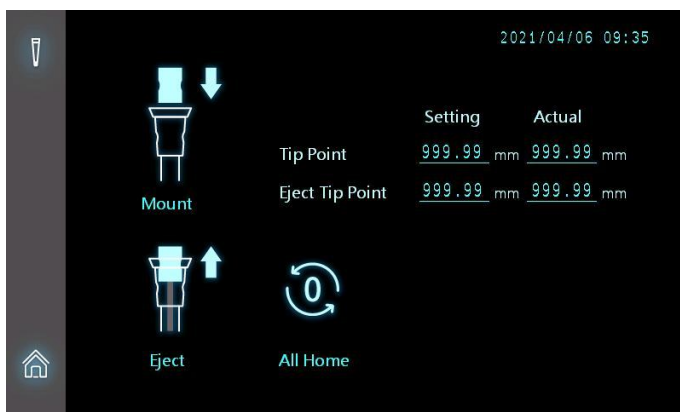- ❑ Tip Point (Setting): Parameter setting value upon mounting tips
- ❑ Tip Point (Actual): Actual tip-mounting parameter value
- ❑ Eject Tip Point (Setting): Parameter setting value upon ejecting tips
- ❑ Eject Tip Point (Actual): Actual tip-ejection parameter


=Remark=

Actual value needs to equal to Setting value

4. User Interface Guidelines

How to test mounting and ejecting spin tips

Step 1: Click "All Home" to set X/Y/Z-axis to zero

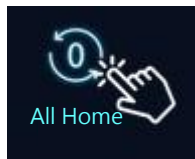Step 2: Place tip combo pack (for use of mounting spin tips) or blank 96DW plate (for use of ejecting spin tips) to #8 plate position


◀ Tip function is designed to test the performance of mounting and ejecting spin tips.

Step 3: Use "▶" and "◀" button to move #8 plate underneath the extraction module

Step 4: Click "Mount" icon to execute mounting spin tips (tip combo pack required) or click "Eject" icon to execute ejecting spin tips (blank 96DW plate required)

Step 5: Check if all spin tips are mounted to spin shaft without evident gap or if all spin tips are completely ejected from spin shaft

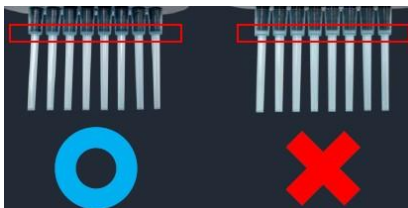
4. User Interface Guidelines

Run (Execution of Extraction Program)


Run function is designed to execute nucleic acid extraction program stored in the instrument.

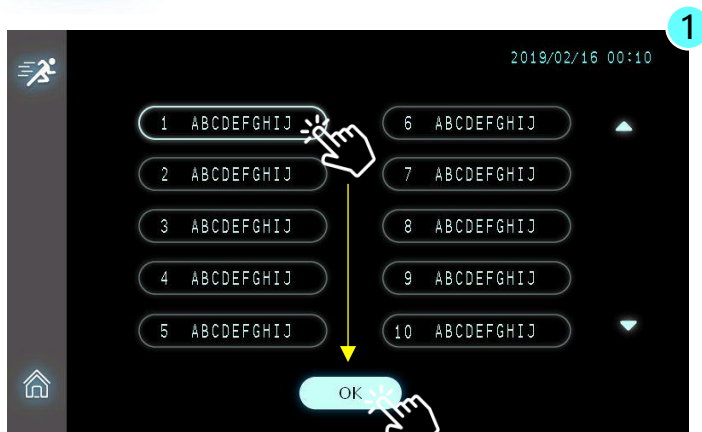- ▲ Maelstrom 9610 is capable of storing up to 25 programs. Please follow the guideline of reagent kits and correctly select the program.


Step 1: Click the program first, and then click [OK] to continue

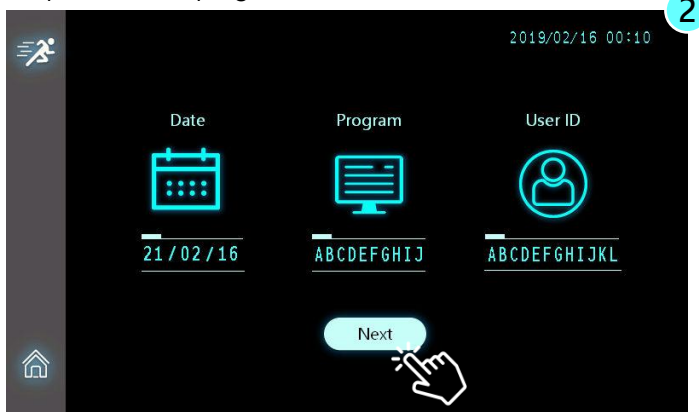
Step 2: Click [Next] to continue

4. User Interface Guidelines

Run (Execution of Extraction Program)

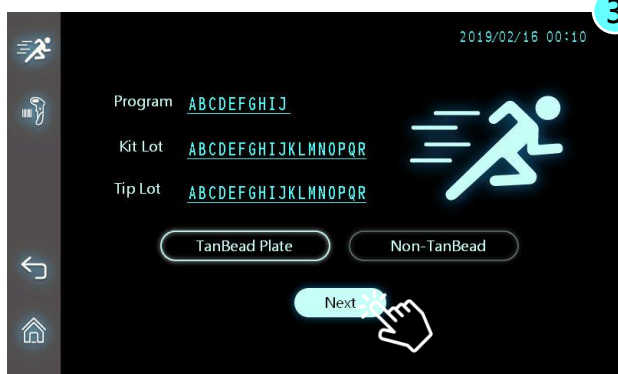Step 3: Click space to input "Kit Lot" and "Tip Lot", if need. Click [Next] to continue. Both Kit Lot and Tip Lot can bear 18 characters input.

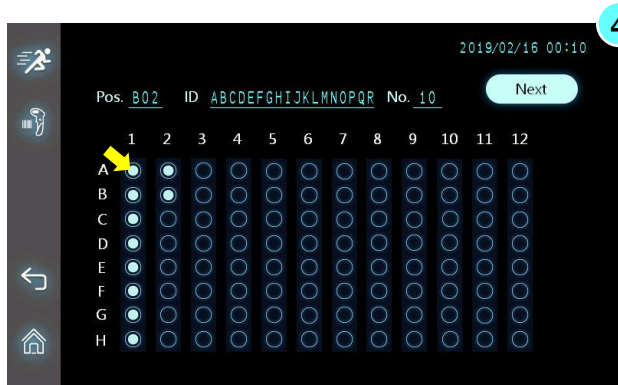Step 4: Input sample ID, if need. Press the circle first, and then use on-screen keyboard or USB barcode scanner to scan in if 1-D barcode is available. Sample ID can bear 18 characters input. Once completed, Press [Next] to continue.


- Pos. : Hole position on the plate
- No. : Total number of Sample ID
- ID : Name of Sample ID

=Remark=

Maelstrom 9610 can be compatible only with USB-type infrared barcode scanner.

4. User Interface Guidelines

Run (Execution of Extraction Program)

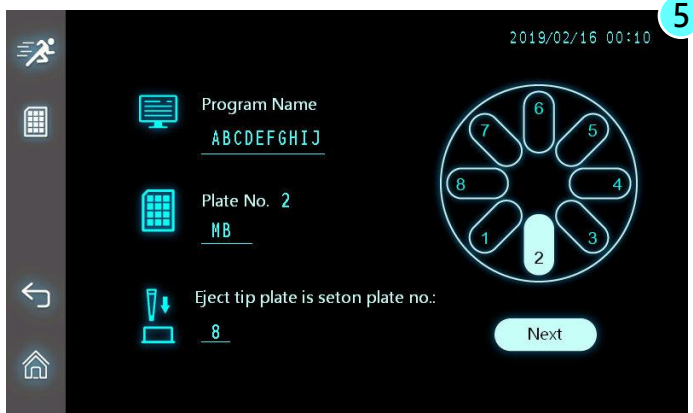Step 5: Follow the guide on screen to place the reagent plate onto correct plate accordingly. Please note that the priority well plate must be spin tips compo pack, please place the eighth plate.

=Remark=

Eject tip plate is allowed to selected from NO.2 to NO.5 plate according to user's needs.


Step 6: Ready to run the extraction program. Click [Yes] to begin the operation. If you choose [No], the process will go back to previous page.

4. User Interface Guidelines

During the Run

Understanding every symbol in the operation process is important. It helps you clearly know the status of nucleic acid extraction.

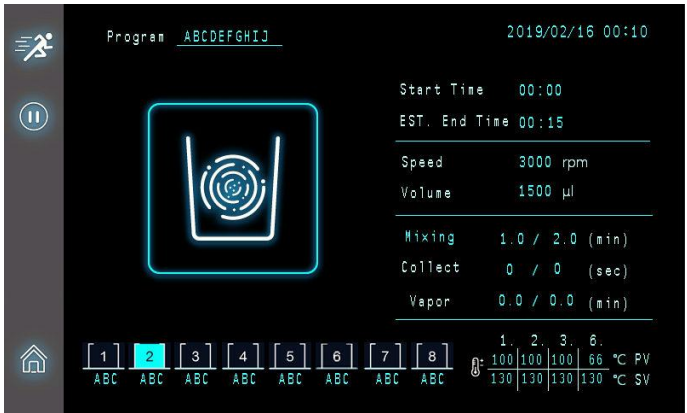Pause: Cease an extraction program temporarily

▲ From plate 1 to plate 6, name of buffer varies from kits to kits

- ❑ Start time: Time of extraction program started
- ❑ Est. End time: Estimated time the extraction program coming to the end
- ❑ Speed: Spin speed
- ❑ Volume: Liquid volume in the working well

4. User Interface Guidelines

During the Run

Understanding every symbol in the operation process is important. It helps you clearly know the status of nucleic acid extraction.

Mixing

Beads collection

In the course of beads collection, whether you click Pause button or Stop button, the instrument will not act immediately until the collection process comes to the end.

Vaporing EtOH

Before entering Elute process, the protocol will must work on Vapor process to evaporate EtOH. Poor evaporation on EtOH will influence the qPCR data significantly.

Program finished

When a run comes to the end, the Buzzer will ring up few seconds.

Resume a run


During the operation, if you click Pause to temporarily stop the program, you can click Resume to continue the operation.

4. User Interface Guidelines


Report (Export extraction records)

Report function enables you to view experiment records and the identity that performed the experiment.

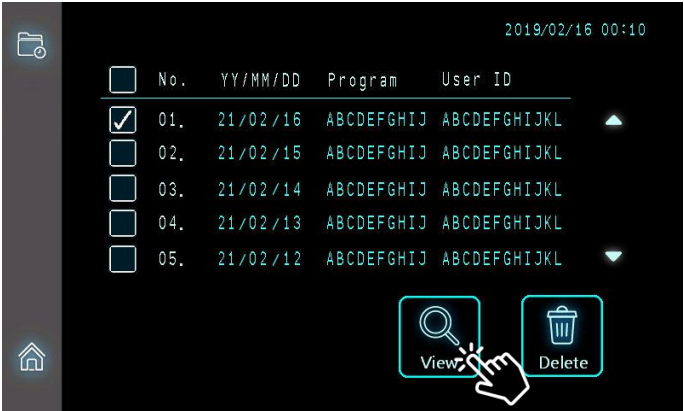▲ Maelstrom 9610 is designed to store the last 15 experiment records.

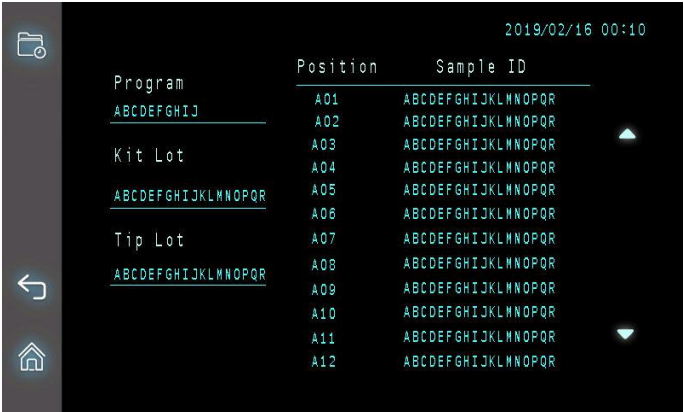
▲ By clicking View button, you can navigate the detail of sample ID.

4. User Interface Guidelines

UV / HEPA

It is highly advised to use UV to clean the instrument after each experiment.

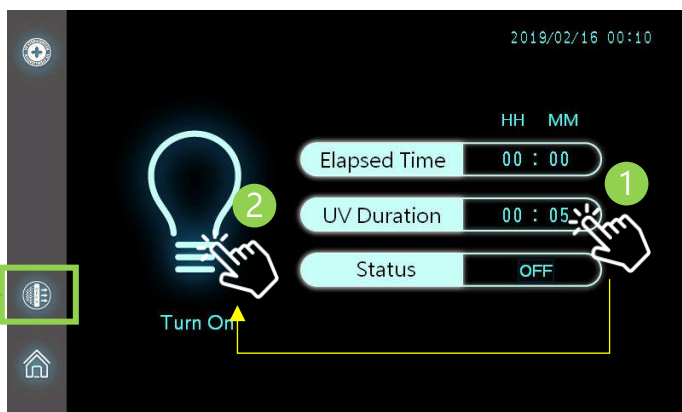▲ Set UV lighting time first, and then press icon to start

- ❑ Step 1: Set lighting time by clicking UV duration field
- ❑ Step 2: Click bulb icon to begin UV lighting

Click HEPA icon to switch to HEPA function
Same with UV, set HEPA time first and then click icon to start


Before turn on UV, please must close the cover door of turntable. The recommended lighting time is at least 30min. The more, the better.

4. User Interface Guidelines

Program (Edit/Delete/Add a Program)

Program function enables you to customize the extraction protocol.

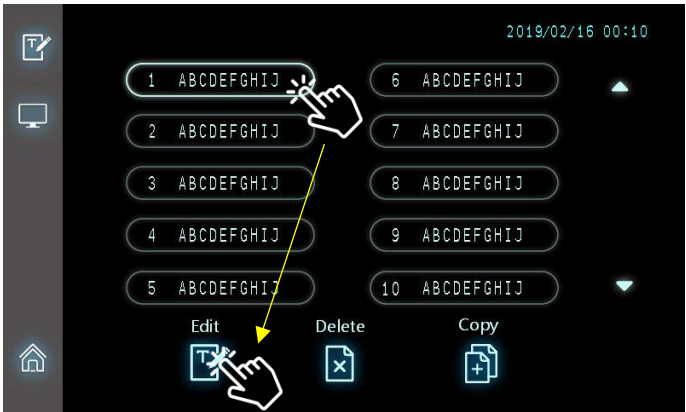▲ Click a protocol first, and then choose a function to proceed

	Edit can be used to adjust the setting of existing protocols, which contains well volume, preheat temp, spin direction, extraction steps, time of mixing and collect.
	To add a new program, you need to use Copy to duplicate a program first, and then use Edit to write the protocol.
	Delete a program no longer used to save the memory space. Maelstrom 9610 can store up to 25 programs and each program can support totally 35 steps.

4. User Interface Guidelines

Edit a Program

Plate	1	2	3	4	5	6	7	8
Volume (μl)	1600	1600	1600	1600	1600	1600	1600	1600
Preheat	100	100	100			100		
Action	----	----	For./Rev.	For.	For./Rev. U/D	For.U/D	----	
Name	ABCD	ABCD	ABCD	ABCD	ABCD	ABCD	ABCD	ABC

Next

▲ Press existing value to change the setting

Program	Name of program, up to 10 characters
Plate	No. of plate
Vol.(μl)	Total processing volume
Preheat (°C)	Heating a plate to a particular temp before mixing step Preheat temp range from 25°C to 100°C
Action	Way of binding during the mixing process For.: Forward For./Rev.: Forward and reverse (spin direction will change regularly) For. U/D: Forward with up-and-down move For./Rev. U/D: Forward and reverse with up-and-down move
Name	Name of plate, which following buffer may be involved in... LB: Lysis buffer WB1: Wash buffer 1 WB2: Wash buffer 2 MB: Magnetic beads EB: Elute buffer

Please note that name of plate 8 is unalterable.

4. User Interface Guidelines

Edit a Program

Step	Extraction step Maelstrom 9610 supports up to 35 steps in a program
Plate	No. of plate
Temp	Please adjust to appropriate protocol temperature before running the test which is suitable for different reagents. Protocol temp and Preheat temp are different. When executing a program, preheat temp will be fulfilled in advance, and protocol temp will be acted as the target plate is about to be used.
Mixing	Time of mixing, range from 0 to 999.9 mins
Spin	Speed of spin (mixing), range from 500 to 3000 rpm
Collect	Time of collecting beads, range from 0 to 999.9 sec
Vapor	Time of evaporate EtOH, range from 0 to 999.9 mins
Pause	Will automatically pause after executing the current step. If additional reagents need to be added during the extraction process, this step can be taken.

Please note that name of plate 1 and plate 8 is unalterable.

4. User Interface Guidelines

Add a new program

1. Touch one of standard programs button.
2. Press Copy button and fill in blank space.
3. Choose new program and press Edit button to set the relevant parameters.

- ▲ Touch one of standard programs button and press Copy button.

- ▲ Confirm the number to add the program. The keyboard will appear after press on destination number. Once complete, click Yes to add the program.

4. User Interface Guidelines

Delete a Program

To delete a program, press the program no longer need and then click Delete icon to execute.

▲ Press the unwanted program, and press the Delete icon


▲ Click Yes to confirm the deletion

4. User Interface Guidelines

Setting

Setting is designed for user to change basic system parameter setting, such as display language, time info and barcode scanner management.

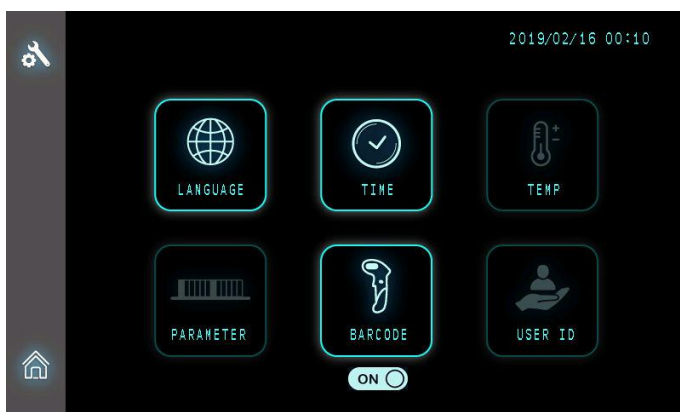- ▲ Admin user is allowed to assess Language, Time, and Barcode functions.

Language

Maelstrom 9610 supports multi-language display, including English, Traditional Chinese, Simple Chinese, Korean and Spanish.

Time

When a new instrument is delivered, please adjust to the right time zone.


Barcode

Scan code to input related information. e.g. lot number, sample ID.

4. User Interface Guidelines

Setting (Language)

Maelstrom 9610 supports multi-language operating interface including English, Traditional Chinese, Simple Chinese, Korean and Spanish.

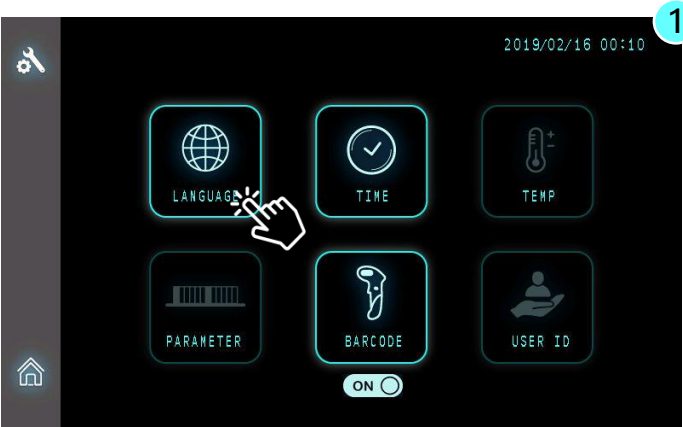▲ Press Language icon to access


▲ Choose the different languages

4. User Interface Guidelines

Setting (Time)

Set correct time for the system according to the local time

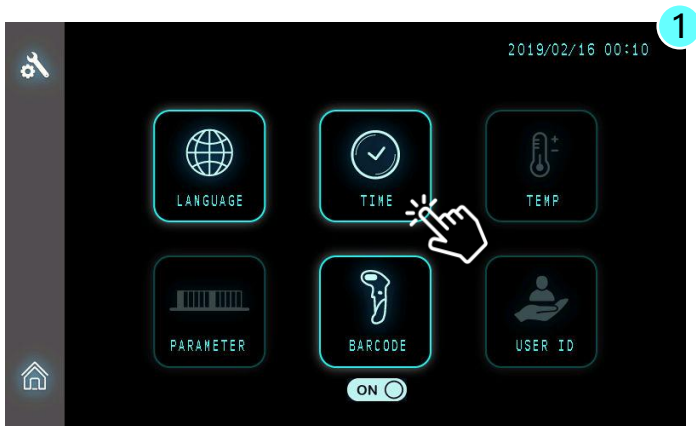
- ▲ Click icon to set Date and Time for the system; once set up is done, press [OK] to apply the latest setting.

4. User Interface Guidelines

Setting (Barcode)

Maelstrom 9610 can be compatible with USB-type infrared barcode scanner.

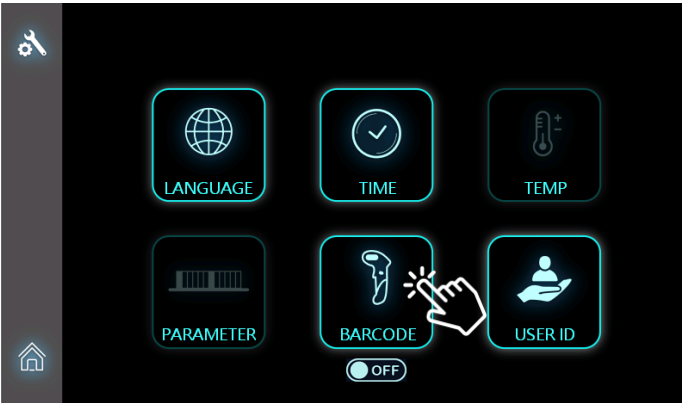
▲ Click the icon to ON/OFF of barcode scanner

⚠ If barcode scanner is set OFF, even USB-type barcode scanner is connected to Maelstrom 9610, scanned data will not be presented.

Setting(USER ID)

The IDs are respectively, Supervisor authority and Operator authority.

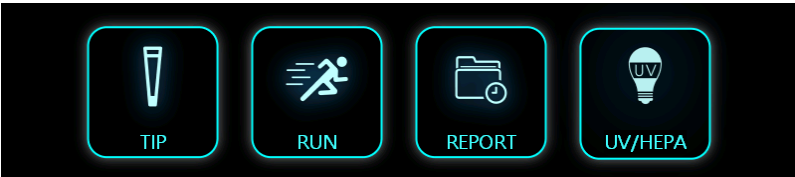▲ Operator authority(four icons)

▲ Supervisor authority(six icons).For Supervisor ID, please contact TANBead.

5. Technical Support

Taiwan Advance Nanotech Inc. provide post-sales service and technical support. In case of any questions, please try to contact our authorized distributor nearest to you or our company for more help.

If there is a need to contact us, you can call our rep number at +886-3-3167568 or email us at service@tanbead.com for further assistance.

To efficiently solve the problems and answer your questions, please provide us your instrument serial number when you talk to our technician personnel. Thanks for your cooperation.

6. Cleaning and Maintenance

- It is important to clean the device after every use. If samples or reagents have been spilled, it is important to clean the device immediately to avoid damage or contamination of samples.
- Wear gloves and appropriate personal protective equipment. If the device is used with biohazardous materials, dispose of any cleaning materials used in accordance with your institutional guidelines.
- The device may go through a run with the magnetic rods unprotected. If this happens, the magnetic rod needs to be cleaned immediately.
- To clean the magnetic rods, wipe with a soft cloth dampened with pure water. Do not use alcohol solvent.
- If the magnetic rods cannot be cleaned, please contact TANBead (Tel: +886-3-3167568) for technical assistance.

7. Disposal

In case the product is to be disposed of, the relevant legal regulations are to be observed. As disposal regulations may differ from country to country, please contact your supplier if necessary. Please follow your institutional and country-specific requirements to handle the disposal of accessories. Device must be decontaminated prior to disposal.

8. Patent


Patent List	
USA	US09616398B2
EU	EP2937136
Canada	CA2862946
Japan	JP6151735B2
Korea	KR101696517B1
China	CN104971638B
Taiwan	TWI526245B
WIPO	WO2016127292

9. About Manufacturer

Here comes the detail of instrument manufacturer and its EC Rep.

- **Manufacturer** : Taiwan Advanced Nanotech Inc.
- **Head Quarter Address**: 6F, No. 188, Wenhe Rd., Guishan Dist., Taoyuan City 333, Taiwan
- **Manufacturer Address** : 4F, No. 188, Wenhe Rd., Guishan Dist., Taoyuan City 333, Taiwan
- **Manufacturer Tel** : +886-3-3167568

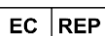- **Contact details** : mdi Europa GmbH
- **Address** : Langenhagener Street 71, D-30855 Langenhagen, Germany
- **Telephone** : +49-511-3908-9530

10. Appendix

Fuse Replacement on Maelstrom 9610

Fuse Specification :

- UL recognized component
- Rated 5A, 250VAC x2pcs
- UFE, Size 5.2Φ x 20mm


Safety precautions of replacing fuses

The following safety precautions will prevent injury to personnel and damage to equipment. These are the MINIMUM safety precautions for replacing fuse.

1. Be sure the power is off in the circuit and the circuit is discharged before replacing a fuse.
2. Use an identical replacement fuse.
3. Remove any corrosion from the fuseholder before replacing the fuse.
4. Be certain the fuse properly fits the fuseholder.

How to fix the fuse

1. Turn the power off and discharge the circuit before removing a blown fuse.
2. Use a fusepuller (an insulated tool) when you remove a fuse from a clip-type fuseholder.
3. Remember, always replace the faulty fuse with one equal in size, amperage, rating and type.
4. Turn the main power switch back on and see if electricity is restored.

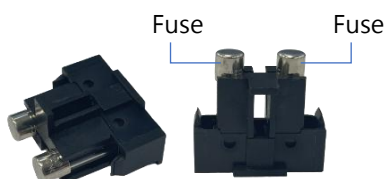
10. Appendix

How to use Maelstrom Master to update extraction protocols in the instrument?


Connect the instrument to your PC

Power on the instrument and then use the cable (USB A to B) to link the instrument to your PC (Type A to PC, Type B to instrument).

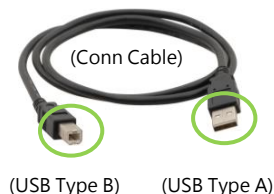
Execute the utility program "Maelstrom 9610 Master"

You can find the utility program from the CD that comes with the instrument. Load the CD to the PC and use Windows Explorer to find the folder named "Maelstrom Master". Then, double click the utility "TANBead Maelstrom 9610 Master".

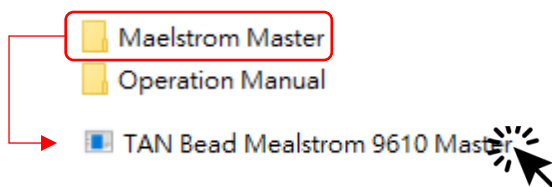The utility program will automatically detect the connection between the instrument and PC. If successfully connected, the pop-out window will tell.

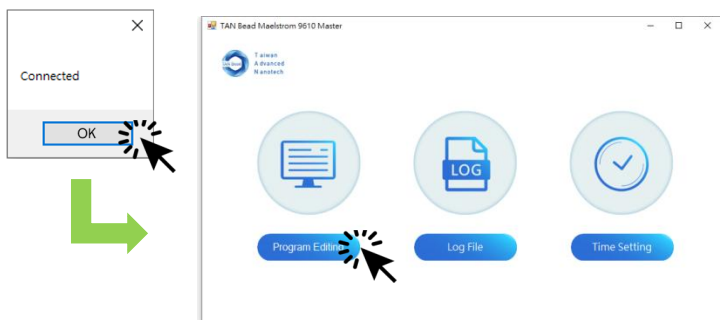
10. Appendix


How to use Maelstrom Master to update extraction protocols in the instrument?

Program Editing

After connection, the utility program will automatically read all protocols stored in the instrument.

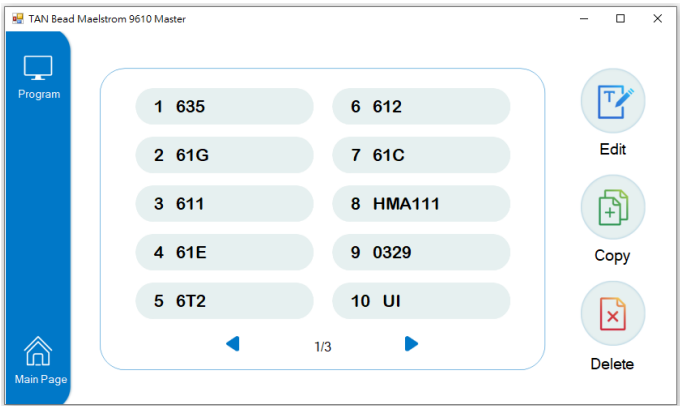Select the protocol first and choose the function to performed. Please note that the protocol you select will be marked in Blue.

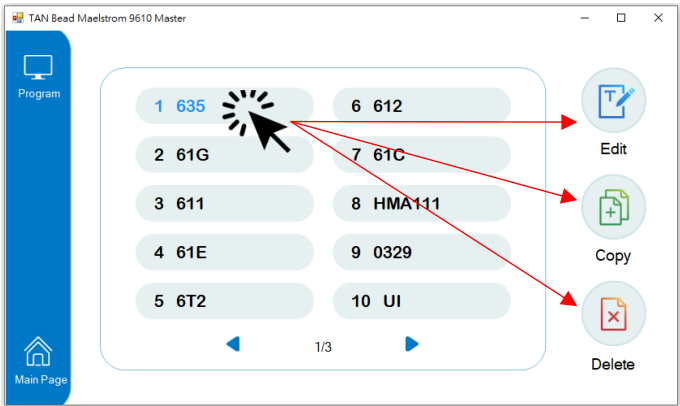
10. Appendix

How to use Maelstrom Master to update extraction protocols in the instrument?

Remember to SAVE settings

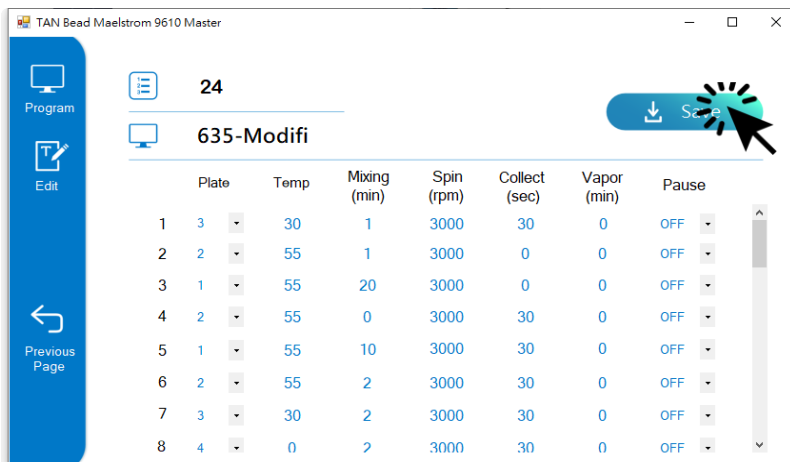After adjustment of protocol contents, please MUST click “Save” to get most updated contents effective. Once you click “Save”, the utility program will automatically update the protocol contents to the instrument.

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